



Win11 Tablet Wifi 16M CMOS Hdmi Microscope Camera

Our Product Introduction

for more products please visit us on cnoec.com

Basic Information

- Place of Origin: China
- Brand Name: CNOEC, OPTO-EDU
- Certification: CE, Rohs
- Model Number: A59.3551
- Minimum Order Quantity: 1 pc
- Price: FOB \$1~1000, Depend on Order Quantity
- Packaging Details: Carton Packing, For Export Transportation
- Delivery Time: 5~20 Days
- Payment Terms: T/T, West Union, Paypal
- Supply Ability: 5000 pcs/ Month




Product Specification

- Sensor: 1/2.3" Sony CMOS Sensor
- Resolution: 16M CMOS
- LCD Screen: 11.6" LCD, Resolution 1920x1080p IPS Screen
- CPU: Inter J4105
- System: Win10 64Bits
- WIFI: 2.4G/5G WIFI, 802.11a/b/g/n/ac
- Output: USB3.0, HDMI, WIFI, Ethernet, Type-C Power Supply
- Highlight: **16M CMOS Hdmi Microscope Camera, Wifi Hdmi Microscope Camera, Win10 Hdmi Microscope Camera**




More Images

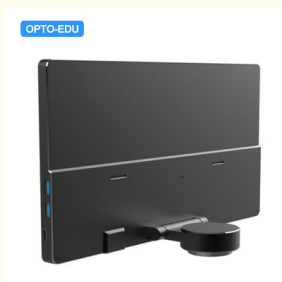



Product Description

16M, 11.6" IPS LCD, 1/2.3" CMOS, LCD Resolution 1920*1080, Win11 64 Tablet, 2.4G/5G WIFI, Bluetooth, USB3.0, HDMI, Ethernet, Type C Power Supply DC12V/2A, Incl. S-EYE Measure Software

Product Details

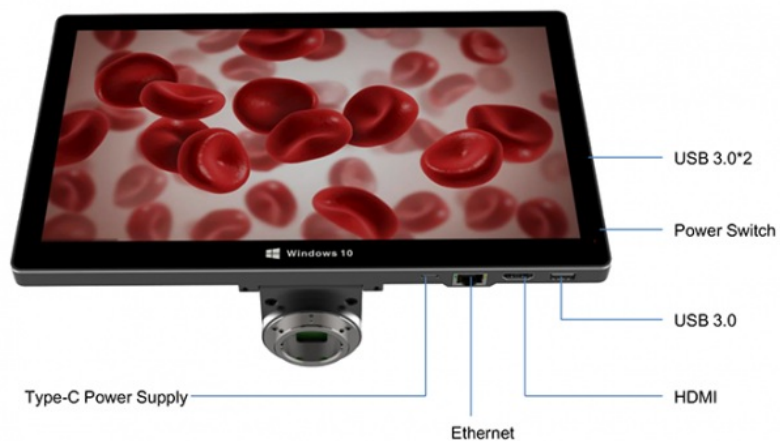

S-EYE Measure Software For Windows

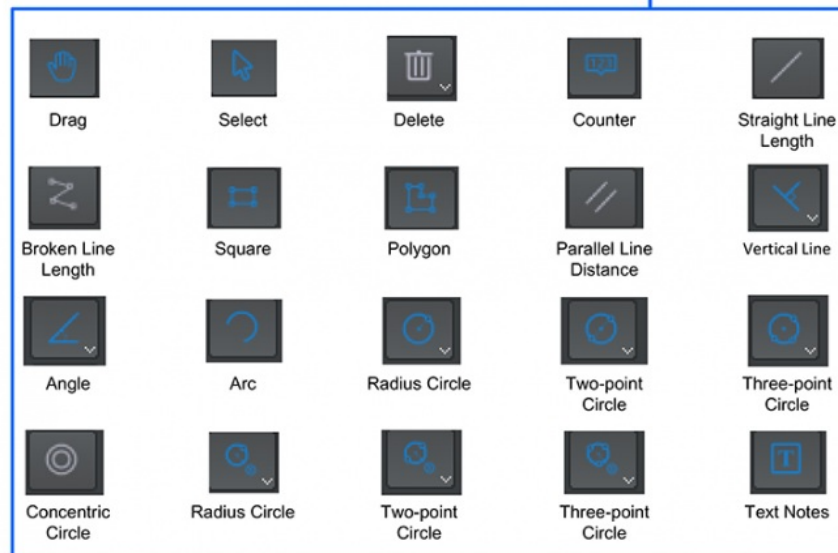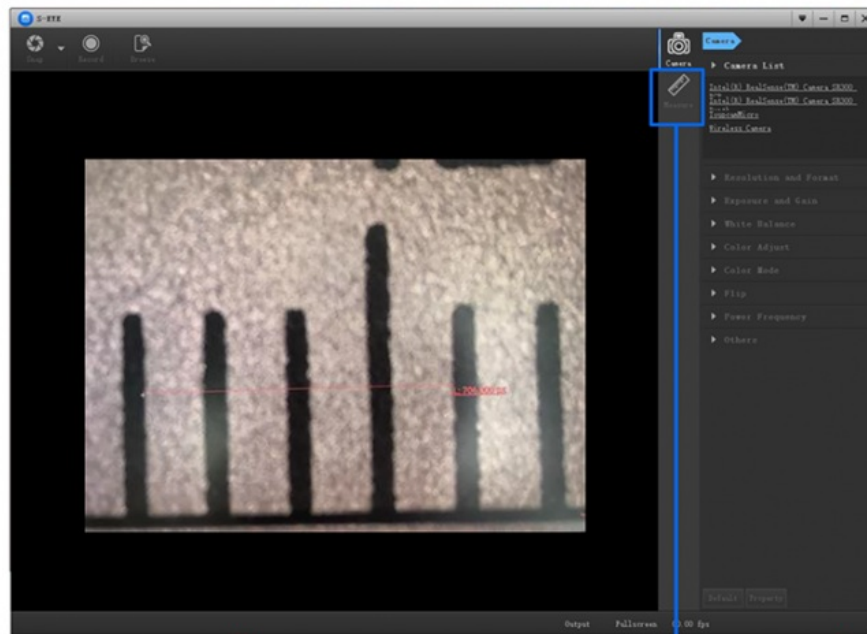

Items	A59.3551 11.6" Win11 Tablet Digital Camera, 5.0M
Sensor	1/2.3" Sony CMOS Sensor
Resolution	16M CMOS
LCD Screen	11.6" LCD, Resolution 1920x1080p IPS Screen
CPU	Inter J4105
System	Win10 64Bits
WIFI	2.4G/5G WIFI, 802.11a/b/g/n/ac
Bluetooth	Bluetooth 4.2
Output	USB3.0, HDMI, WIFI, Ethernet, Type-C Power Supply
Pixel Size	1.335x1.335um
Power Supply	DC 12V/2A Power Adapter, Input 100-240V AC
Function	Photo 16M, Video 1080p Full HD
	Adjust: Brightness, Contrast, Saturation, White Balance, Exposure, Flicker Suppression, Sharpness, Gamma, CB/Cr Gain
	UI Interface: Mouse Control, Touch Screen
	Storage Photo, Video, Preview, SD Card Format

Measure	<p>S-EYE Measure Software:</p> <p>Length: Straight Line Length, Broken Line Length, Curve Length, Parallel Line Distance, Point Line Distance, Line Length</p> <p>Geometry: Radius Circle, Two-point Circle, Three-point Circle, Concentric Circle, Radius Circle, Two-point Circle, Three-point Circle</p> <p>Geometric Area: Polygon, Square</p> <p>Remark: Point Coordinates, Crosshairs, Coordinate Systems, Text Notes</p>
Size / Weight	Product Size 285x177x12.5mm, Net Weight 0.9kg
Standard Package	LCD Camera, DC12V Power Adapter
Quick Manual	<p>1. First, connect the product to Type-C DC-12V power supply</p> <p>2. After connecting the power supply, press the power on button, and the blue light of the product will always be on.</p> <p>3. Double click the S-EYE Measure Software on the desktop of Window 11</p> <p>4. Click "Open Video Source" in S-EYE software, and choose camera YW1600</p> <p>6. Measure and adjust as you need.</p>



Opto-Edu (Beijing) Co., Ltd.



0086 13911110627



sale@optoedu.com



cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China