for more products please visit us on cnoec.com

OPTO-EDU 3.2x - 192x Binocular Forensic Comparison Microscope with LED Illumination for Student and School Applications

Basic Information

Place of Origin: China

Brand Name: CNOEC, OPTO-EDU

Certification: CE, Rohs
Model Number: A18.1825
Minimum Order Quantity: 1 pc

Price: FOB \$1~1000, Depend on Order Quantity
 Packaging Details: Carton Packing, For Export Transportation

• Delivery Time: 5~20 Days

Payment Terms: T/T,West Union,Paypal
 Supply Ability: 5000 pcs/ Month



Product Specification

Applications: Student|school
 Certification: CE|Rohs
 Drawtube: Binocular
 Illumination: LED
 Magnification: 3.2x - 192x
 Theory: Comparison Mic

Theory: Comparison Microscope
 Highlight: binocular optical microscope, bullet comparison microscope



More Images











Product Description

OPTO-EDU 3.2x - 192x Binocular LED Bullet Forensic Comparison Microscope A18.1825

Forensic Comparison Microscope

The A18.1825 comparison microscope, newly developed in 2013, is designed for forensic ballistics applications and education. This instrument provides students and professionals with a complete comparison microscope system at a competitive price point, featuring steady quality and full functionality.

Key Features

- Total Magnification: 3.8x~230x
- Optical System: Wide Field Eyepiece WF10x/22mm, WF20x/12mm with High Eyepoint
- Zoom Objective: 5 Step Zoom Plan Objective (0.8x, 1.25x, 2x, 3.2x, 4.8x)
- Auxiliary Lenses: Includes 0.4x and 2.0x lenses
- Comparison Modes: Full left/right, split-image or overlapping image comparison
- Advanced Systems: Magnification Correcting System with Correcting Lens
- Adjustment: Separation Line Adjusting System with special metal adjust driver
- Dual Lighting: 12V50W fan-cooled halogen light and 5W LED cold light, both adjustable

- Control Panel: Touch panel switch for separate or combined light source control
- Comprehensive Accessories: Includes coaxial light, polarizing light, UV lens, and color filters
- Power: 80~240V wide range voltage



A18.1825 Digital Forensic Comparison Microscope




Product Details



Large Binocular Head both eyepiece tube diopter+/-5 adjustable, wit WF 10X/22mm eyepiece high eyepoint




Magnification Switch
5 step zoom plan objective: 0.8x, 1.25x, 2x, 3.2x
4.8x

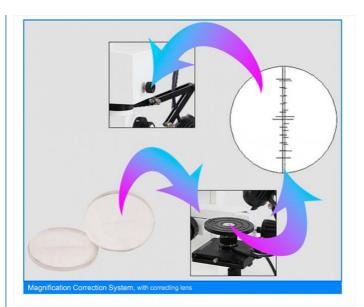







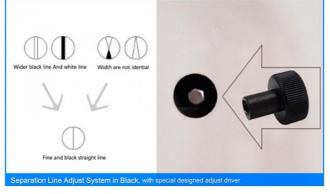









Aluminum Box Package



Product Specification

Eyepiece Mag./Viewfield		Eyepiece10X	Eyepiece20X	Photo Adapter	Video Adapter	Working
Objective	Bridge Body					
0.8X	1.2X	9.6X/φ28	19.2X/φ17	2.5X	0.5X	101mm
1.25X		15X/φ18	19.2X/φ17			
2X		24Χ/φ11	48X/φ7			
3.2X		36Χ/φ7	77X/φ4.5			
4.8X		58X/φ4.2	115Χ/φ2.3			

Technical Specifications

Specification	Details
Magnification	3.8x~230x
Head	Binocular Head, 30° Incline, Interpupilary Distance 48~75mm, Diopter +/-5 Adjustable
Eyepiece	WF10x/22mm, High Eyepoint WF20x/12mm, High Eyepoint
Objective	5 Step Zoom Plan Objective 0.8x, 1.25x, 2x, 3.2x, 4.8x Auxillary Lens 0.4x, 2.0x
Working Distance	101mm
View Field	Dia. 2.3~62mm
Magnification Correcting	Magnification Correction System, With Correcting Lens
Separation Line Adjusting	Separation Line Adjust System In Back, With Special Designed Adjust Driver
Working Stage	Square Manual Working Stage 90°90mm, Moving Range 54°54°54″54mm Round Titable Working Stage Dia. 90mm, 0-45° Inclined, 360° Rotatable Joint Left/Right Moving Range 54mm, Joint Up/Down Range 60mm
Light Source	12V50W Fan Cooling Strong Halogen Light (2pcs) 5W LED Cold Light (2pcs) Coaxial Light Lens (2pcs) Polarizing Light Lens (2pcs) Ultraviolet Light Lens (2pcs) Filters: Red, Blue, Yellow, Green (2pcs each)
Adapter	0.5x CCD Adapter For Digital Camera (1pc) Front View Photo Adapter 2.5X For SLR Digital Camera (Nikon And Canon)
Power	80~240V Wide Range Voltage
Bullet Holder	Standard Bullet Holder 2 Sets
Dimension	Instrument Size 37*30*63cm

Optional Accessories

Accessory	Details
Light Source	150W Cold Light Source, With Optical Fiber & Light Collector (A56.1821)
Bullet Holder	Motorized Bullet Holder 1 Set (A54.1813)
LCD Pad Digital Camera	9.7" LCD Pad Digital Camera, Android System (A59.3520)

More Pictures











Opto-Edu (Beijing) Co., Ltd.

© 0086 13911110627 Sale@optoedu.com cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China