



## 1000x Metal Polarizing Light Microscope Halogen Lamp Microscopes A15.0201

### Our Product Introduction

#### Basic Information

- Place of Origin: China
- Brand Name: CNOEC, OPTO-EDU
- Certification: CE, Rohs
- Model Number: A15.0201
- Minimum Order Quantity: 1 pc
- Price: Negotiation
- Packaging Details: Carton Packing, For Export Transportation
- Delivery Time: 5~20 Days
- Payment Terms: T/T, West Union, Paypal
- Supply Ability: 5000 pcs/ Month



#### Product Specification

- Application: Research|metal
- Certification: CE|Rohs
- Drawtube: Binocular
- Illumination: Halogen Lamp
- Magnification: 40x - 1000x
- Theory: Polarizing Light Microscope
- Name: 1000x Metal Polarizing Light Microscope  
Halogen Lamp Microscopes A15.0201
- Keyword: Polarizing Microscope
- Highlight: 1000x polarizing light microscope ,  
metal polarizing light microscope ,  
halogen lamp polarizing microscope



## Product Description


### Description:

A15.0201 polarizing microscope is utilized polarize light to observe and research the polarize characteristic of the matter which have the double refraction .It can be used to medical inspection, such as the frame detail of the nerve fibre,the animal muscle, the plant fibre and the crystal of the tooth,bone,hair and alive cell etc. It can be observe the fine frame of the abiochemistry salt crystal that the fine frame can't be observe for nature light.

### Specifications:

A15.0201 Polarizing Microscope Specification		AB
Eyepieces	WF10X/18mm Wide Field	●●
Objectives	Achromatic 4X/0.10	●●
	Achromatic 10X/0.25	●●
	Achromatic 40X/0.65 (Spring)	●●
	Achromatic 100X/1.25(Spring,oil)	●●
Head	Binocular inclination of 30°	●●
Analyzer	Sliding	●●
Focus System	Coaxial coarse/fine system, with tension adjustable and up stop, minimum division of fine focusing 2um	●●
Nosepiece	Quadruple, backward ball bearing inner locating	●●
Stage	Rotatable stage, diameter 120mm	●
	Double layer mechanical stage135x125mm, moving range 75x35mm, with rotatable stage dia.100mm	●
Condenser	Abbe Condenser N.A. 1.25 rack & Pinion adjustable	●●
Filter	Frosted filter	●●
	Blue filter	●●
Polarizer	360° rotatable stage	●●
Collector	For illumination with halogen lamp (with field diaphragm)	●●
Illumination Unit	With polarizer, 6V20w, halogen lamp, adjustable brightness	●●
Optional Accessories		Item No.
Eyepiece	WF16x/11mm wide field	A51.0203-16A
	Dividing WF10x/18mm, 0.1mm/Div	A51.0205-10
Head	Monocular, 360°, inclination of 30°	A53.0201-5
	Trinocular, inclination of 30°, observe and photograph at same time	A53.0207-2
Nosepiece	Quintuple, backward ball bearing inner locating	A54.0203-1
Filter	Green filter	A56.0209-2
	Yellow filter	A56.0209-3
CCD Adapter	0.4X	A55.0202-01
	0.5X	A55.0202-04
	1X	A55.0202-02
	0.5X with dividing 0.1mm/Div	A55.0202-03
Photo Adapter	2.5X/4X Change over photograph adapter with 10X viewing eyepiece	A55.0201-1
	4X Focusing photograph adapter	A55.0201-2
	MD Adapter	A55.0201-3
	PK Adapter	A55.0201-4
Adapter	Digital camera adapter CANON (A610,A620,A630,A640)	A55.0204-01



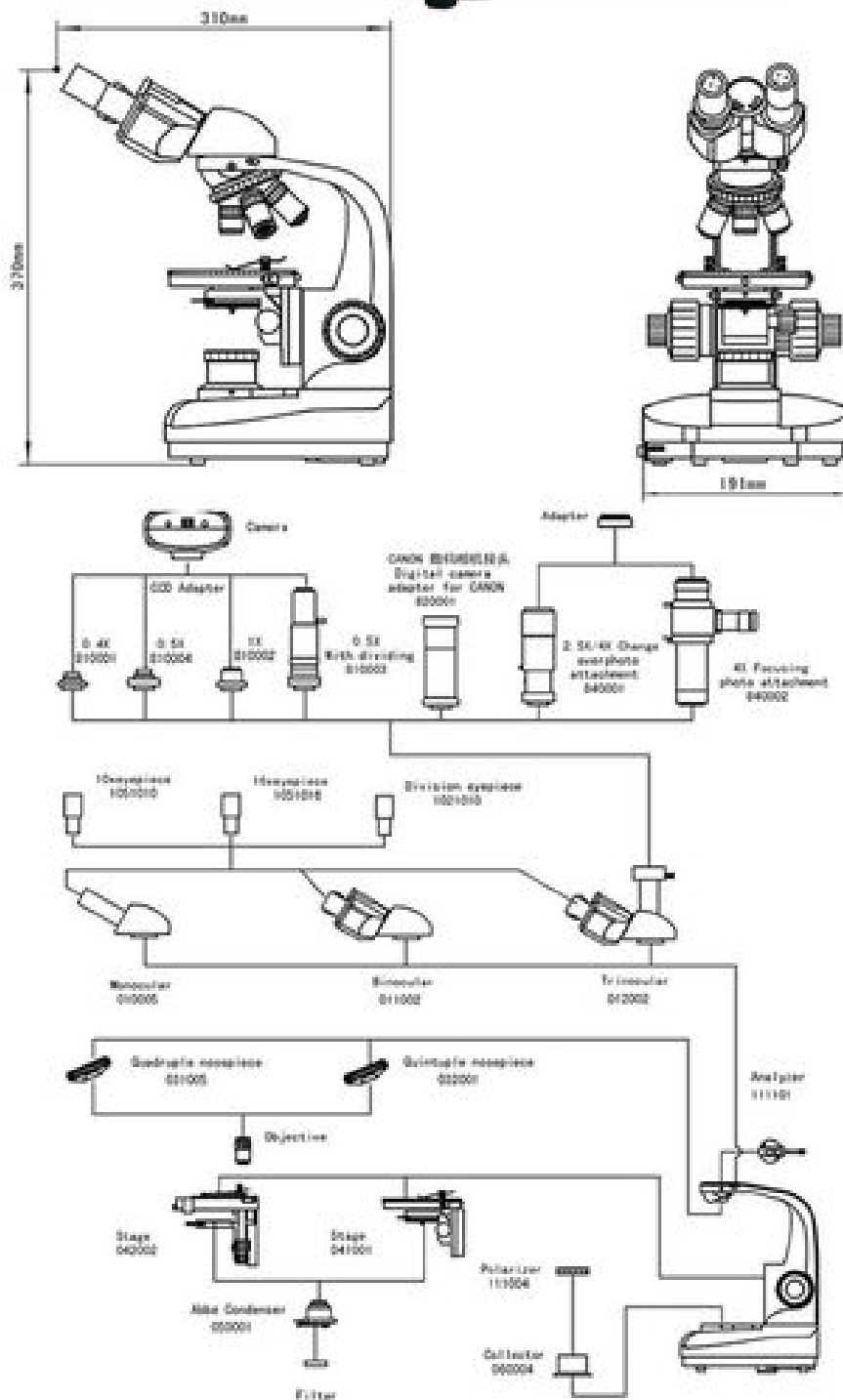

**Opto-Edu (Beijing) Co., Ltd.**



0086 13911110627



sale@optoedu.com



cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China