



OPTO-EDU A13.2605-A Metallurgical Microscope, Binocular, Reflect Light

Our Product Introduction

Basic Information

- Place of Origin: China
- Brand Name: CNOEC, OPTO-EDU
- Certification: CE, Rohs
- Model Number: A13.2605-A
- Minimum Order Quantity: 1 pc
- Price: FOB \$1~1000, Depend on Order Quantity
- Packaging Details: Carton Packing, For Export Transportation
- Delivery Time: 5~20 Days
- Payment Terms: T/T, West Union, Paypal
- Supply Ability: 5000 pcs/ Month



Product Specification

- Certification: CE
- Drawtube: Binocular|Trinocular
- Illumination: LED
- Highlight: **usb digital microscope,**
infinity plan microscope



Product Description

Metallurgical Microscope, Binocular, Reflect Light

Entry Level Economical High Quality Upright Metallurgical Microscope
 Wide Field Eyepiece WF10x/20mm + Reticle Eyepiece With Micrometer 0.1mm
 Metallurgical LWD Infinity Plan Objective 5x,10x,20x,50x,80x,100x(Dry) Optional
 Mechanical Working Stage 200x150mm, Moving Range 76x52mm With No Hole

Reflect Light Source 5W LED With Five Holes Filter Disc

Items	A13.2605 Metallurgical Microscope Specification	A	B	Cata. No.
Head	Seidentopf Binocular Head Inclined 30°, Eyetube Rotatable 360°, Interpupillary Distance 48-75mm	●		A53.2621-B
	Seidentopf Trinocular Head Inclined 30°, Eyetube Rotatable 360°, Interpupillary Distance 48-75mm, Light Split Fixed P20/E80		●	A53.2621-T
Eyepiece	WF10x/20mm, High Eyepoint,	●	●	A51.2602-1018
	WF10x/20mm, With Micrometer X-0.1mm*180 Div / Y-Single Line, Diopter Adjustable	●	●	A51.2605-1018D
Nosepiece	Quadruple, Backward	●●	●●	A54.2610-N734
Metallurgical LWD Infinity Plan Objectives	LPL 5x/0.13, W.D.=16.04mm	●	●	A5M.2633-5
	LPL10x/0.25, W.D.=18.48mm	●	●	A5M.263310
	LPL 20x/0.40, W.D.=8.35mm	○	○	A5M.2633-20
	LPL 50x/0.70, W.D.=1.95mm	●	●	A5M.2633-50
	LPL 80x/0.80, W.D.=0.85mm	○	○	A5M.2633-80
	LPL 100x/0.90(Dry), W.D.=1.00mm	○	○	A5M.2633-100
Working Stage	Double Layer Mechanical Stage Size 200×140mm, Travel Range 76×52mm	●	●	A54.2601-S04T
Focusing	Coaxial Coarse And Fine Focusing Knobs, Travel Range: 22mm, Scale: 2μm	●	●	
Illumination	Reflect Light Source, 5WLED Warm Lamp, With Five Holes Filter Disc (Blue, Green, Amber, Grey, Empty) Powered By Safety Low Voltage Charger, Input Wide Voltage 100V-240V, Output 5V1A, Support Power Bank Supply For Outdoor Use	●	●	A56.2613-6V30W
Polarizing	Analyzer + Polarizer	○	○	A5P.2603-BS
Adapter	C-Mount 0.5x, Focus Adjustable	○	○	A55.2601-05
	C-Mount 0.75x, Focus Adjustable	○	○	A55.2601-75
	C-Mount1.0x, Focus Adjustable	○	○	A55.2601-10
	Digital Eyepiece Camera Adapter, Dia. 23.2mm/25mm	○	○	A55.2610
Package	Carton Size 570*275*450mm, 1 Unit/Ctn			
Weight	Gross Weight 10.5kgs, Net Weight 8.0kgs			
Note: "●" In Table Is Standard outfits, "○" Is Optional Accessories.				



OPTO-EDU (BEIJING) CO., LTD.

F-1501 Wanda Plaza, No.18 Shijingshan Road, Beijing 100043, China
 Tel:+8610 88696020 Fax:+8610 88696085

A13.2605

Metallurgical Microscope, Reflect Light




Product Models & Size(mm)



A13.2605-A



A13.2605-B

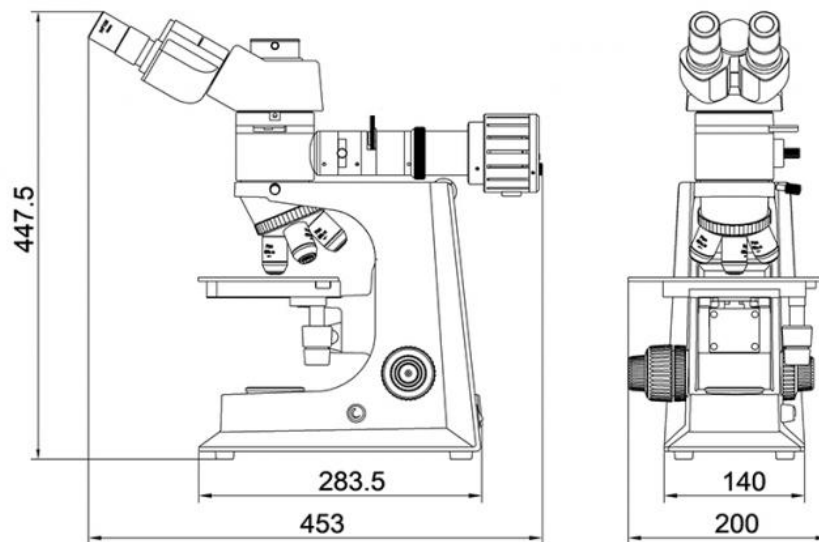


A13.2605 Inverted Metallurgical Microscope Size(mm)

Product Conformation & Accessories

- 01-Eyepiece
- 02-Binocular (Trinocular) Observing Tube
- 03-Locking screw for Camera
- 04-Locking screw for Observing Tube
- 05-Filter turret
- 06-Aperture diaphragm adjusting ring
- 07-Light Box
- 08-Polarizer
- 09-Locking screw for main body
- 10-Objective
- 11-Right coarse focusing hand wheel
- 12-Right fine focusing hand wheel
- 13-Dimming potentiometer
- 14-Main body
- 15-Stage
- 16-Mechanical stage
- 17-Nosepiece
- 18-Epi illumination device
- 19-Analyzer



A5P.2603-BS Analyzer + Polarizer

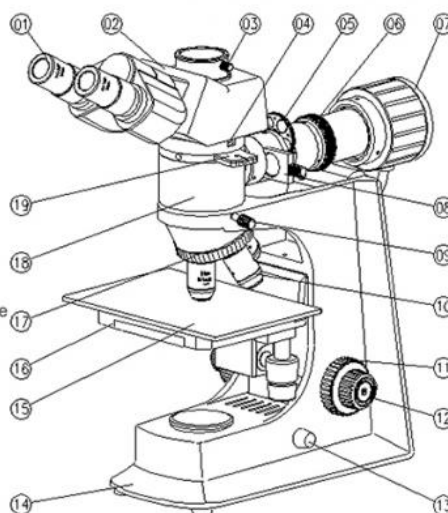

**A55.2610
Digital Eyepiece
Camera Adapter**



**A55.2601-05
C-Mount 0.5x,
Focus Adjustable**



**A55.2601-75
C-Mount 0.75x,
Focus Adjustable**



**A55.2601-10
C-Mount 1.0x,
Focus Adjustable**



**A5M.2633
Metallurgical LWD Infinity Plan Objectives**



Opto-Edu (Beijing) Co., Ltd.

☎ 0086 13911110627

✉ sale@optoedu.com

🌐 cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China