



WIFI 12X-65X Handheld Digital Microscope For iPhone / iPad / PC / Android A34.4902

Our Product Introduction

for more products please visit us on cnoec.com

Basic Information

- Place of Origin: China
- Brand Name: CNOEC, OPTO-EDU
- Certification: CE, Rohs
- Model Number: A34.4902
- Minimum Order Quantity: 1 pc
- Price: FOB \$1~1000, Depend on Order Quantity
- Packaging Details: Carton Packing, For Export Transportation
- Delivery Time: 5~20 Days
- Payment Terms: T/T, West Union, Paypal
- Supply Ability: 5000 pcs/ Month



Product Specification

- Highlight: Portable Digital Microscope, Hand held Digital Microscope



More Images



Product Description

WIFI 12X-65X Handheld Digital Microscope For iPhone / iPad / PC / Android A34.4902

A34.4902 Wireless WIFI USB Digital Microscope, 12~65x, For iPhone/iPad/PC/Android, Newly Designed in 2014, is your best tool to view microworld on your mobile phone, pad or computer! It works on its built-in Wi-Fi module, it allows up to about 200 mobile devices like iPhone, iPad, Android phones, PC Tablet to capture images/video simultaneously. The handheld digital microscope KoPa W5, a portable and light-weighted Wi-Fi microscope that fit in the palm of your hand, providing remarkable image on precision inspection, skin care, education, art preservation, science, you will find dozens of uses for the Wi-Fi microscope W5 with up to 65X magnify.

Watch and share tiny but high-resolution images/videos in real time!

1) Wireless sharing:

Captured images/video can be viewed on multiple mobile device (iPhone/iPad, Android phones, PC with Wi-Fi function, tablets). Supposed there are 5 mobile Wi-Fi devices need to connect with the Wi-Fi microscope, each users could simultaneously see what the microscope captured by download free live camera software in APP. If there are more than one A34.4902 in a room, we suggest changing their names for distinction.

2) USB sharing:

For facilitate your work, the A34.4902 provide you two working modes: Wi-Fi and USB working modes. Image you are inspecting an antique/ art painting/signature/insect, distinguish fake money, you may use plug the A34.4902 with normal USB wire into computer and share inspected content with your colleagues/friends/family on Skype,MSN,QQ,or other instant chatting tools

Application:

The Wi-Fi microscope A34.4902 has widely application, we accept ODM, design customizing hardware and software according to your industrial requirement.

- Health care: for instance, dental care. You may check your dental conditional by using the microscope at any time.
- Industrial application: for example, PC board enlarges checking.
- Valuables discrimination: for example, banknote, jewelry, and antique discrimination.
- Quality control Textile, leather products, oil pipeline inspection etc.
- Indoor wireless observation /Wireless wilderness exploration/ Collections /hobbies/ skin/hair

A34.4902 Wireless WIFI USB Digital Microscope, 12~65x, For iPhone/iPad/Pc/Android	
Magnification	12~65x
Sensor	5.0M True Color CMOS
Interface	USB 2.0, No Need For External Power
Resolution	2592*1944 Max.
Frame Rate	40 Fps, VGA
White Balance	Automatic
Exposure	Automatic
Data Format	Image: BMP, PNG, JPEG, GIF Video: AVI
Wireless Protocol	IEEE 802.11B/G/N
Frequency	2.400~2.483GHz, According to Local Regulation
Transmission	150Mbps Max.
System Support	iOS of iPad, iPhone; Android, Windows 32 & 64 Bits
Illumination	Built-in 8 High Light LEDs
Power	Rechargeable Lithium Battery 1200mA, Self Protection, Working Time 120 Minutes
Dimension	60mm Diameter *54mm Height
Net Weight	100g
Packing	Color Packing Box 8*9*10cm, 18 pcs/Carton, Size 61*24*20cm, G.W.6.5Kgs
4 Ways To Watch	A. Watch the images with Safari /Google Chrome explorer website: http://10.10.1.1:8196 (defalut IP address :10.10.1.1 for iOS 5) B.For iPhone,iPad,iPod touch,iOS Download free APP "Feboncam" , "Feboncam"to view pictures Note:FEBON CAM only support ipad. C. For Android Download free software"IP_Cam View Lite" D.For PC Windows 1)Use our OSTEC own software "KoPa capture" 2)Download free software name "AMCap" or " HiCap"



 A34.4902 Wireless USB Digital Microscope, 12~65x, For iPhone/iPad/Pc/Android



A34.4902 Wireless USB Digital Microscope, 12~65x, With Skin View Stand



A34.4902 Packing Gift Box

Including Microscoe & Skin View Stand






A34.4902 Wireless USB Digital Microscope, 12~65x, Support USB 2.0 Interface



Used With iPad, For Teaching, Exploration, Collections, Hobbies, Skin, Hair Inspection



A34.4902 Different Stands for Different Objects Viewing

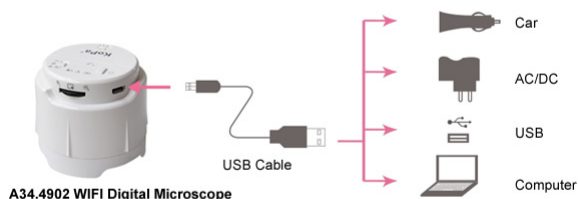

A34.4902 Digital Microscope




1. Turn On Microscope
2. Turn On WIFI Devices
3. View, Capture, And Record

A34.4902 Easy To Use

A34.4902 LED Light Adjustable

Simply Press Button On Top, You Can Adjust LED Light From Off to Full By 8 Levels

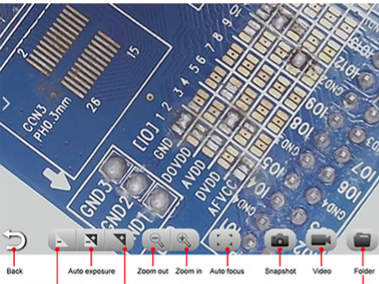

A34.4902 Ucam App For IOS, Android

Ucam is an independent developed App dedicated for Ostec WIFI imaging camera products such as WIFI microscope, WIFI spotting scope, WIFI eyepiece, etc. It has strong functions featured with snapshot, record video, auto focus, folder management, and it is super easy to operate.





Search Ucam App at APP store (IOS users), Google play (Android users), and install the App. Turn on WIFI work mode of Ostec WIFI imaging camera products, when indicator light start flashing, click settings on your IOS/Android devices, and connect with your Ostec WIFI imaging camera.

SSID : VideoWIFI-XXXXXX

Default passcode for connecting : 12345678


Default passcode for opening uCam: ostec

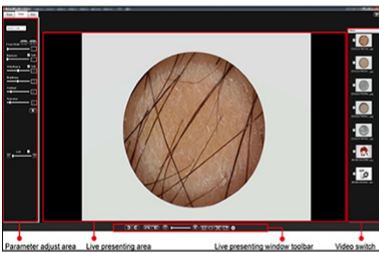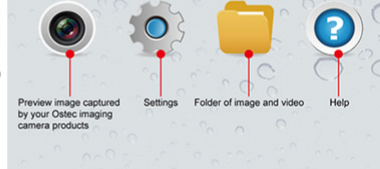

A34.4902 Micro Cam Software For Windows

Supplied in CD-ROM With Microscope It can preview real-time image, capture and record easily.

Micro Cam supply full range adjustment for all parameters, enable you to get the best shape image quality as needed.

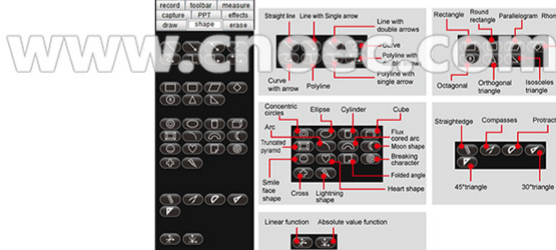



Opto-Edu (Beijing) Co., Ltd.

☎ 0086 13911110627

✉ sale@optoedu.com

🌐 cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China