



Laboratory Research Monocular Polarizing Light Microscope CE A15.1016

Our Product Introduction

Basic Information

- Place of Origin: China
- Brand Name: CNOEC, OPTO-EDU
- Certification: CE, Rohs
- Model Number: A15.1016
- Minimum Order Quantity: 1 pc
- Price: Negotiation
- Packaging Details: Carton Packing, For Export Transportation
- Delivery Time: 5~20 Days
- Payment Terms: T/T, West Union, Paypal
- Supply Ability: 5000 pcs/ Month



Product Specification

- Application: Research, Laboratory
- Certification: CE
- Drawtube: Binocular
- Illumination: Halogen Lamp
- Theory: Polarizing Light Microscope
- Highlight: **student compound microscope ,
compound light microscope**



Product Description

Laboratory Research Monocular Polarizing Light Microscope CE A15.1016

A15.1016 Polarizing Microscope, With features of easy operation, complete functions and economy, it is an ideal instrument in geology, petroleum, coal, chemicals and physics fields. It is also widely used in academic demonstration and research.

- Monocular, Binocular, Trinocular Head, 45° Inclined, 360° Rotatable
- Strain Free Achromatic Objective
- Round Stage Dia. 135mm, 360° Rotatable
- Rotatable Analyzer With Optical Compensator
- Swing-Out Polarizer 0~90° Rotatable, Located On Top Of Collector

A15.1016 Polarizing Microscope		M	B	T
Head	Monocular Head, 45° Inclined, 360° Rotatable	o		
	Sliding Binocular Head, 45° Inclined, 360° Rotatable		o	
	Sliding Trinocular Head, 45° Inclined, 360° Rotatable			o
Eyepiece	WF10x/18mm		o	o
	WF10x/18mm, With Scale Of Crosshair			
Objective	Strain Free Achromatic Objective 4x, 10x, 40x			
Nosepiece	Outward Quadruple Nosepiece			
Analyzer	Rotatable Analyzer With Graduation 0~90° Bertrand Lens, Sliding In/Out Of Optical Path Optical Compensator : -- λ Slip, First Class Red -- $1/4\lambda$ Slip -- Quartz Wedge			
Working Stage	Round Stage Dia. 135mm, 360° Rotatable, Graduation 1° Increment, Minimum Resolution 6', Center Adjustable			
Focusing	Coaxial Coarse & Fine Adjustment, Focusing Range 28mm, Fine Focusing Division 0.002mm			
Condenser	Abbe Condenser N.A.1.25, With Iris Diaphragm & Filter			
Polarizer	Swing-Out Polarizer 0~90° Rotatable, Located On Top Of Collector			
Illumination	Halogen Lamp 6V20W, Brightness Adjustable			



A15.1016 Polarizing Microscope

Rotatable Analyzer With Gradation 0°-90°
-- λ Slip First Class Red
-- $1/4\lambda$ Slip Quartz Wedge

www.cnoec.com

Round Stage $\Phi 135\text{mm}$, 360° Rotatable




-- λ Slip First Class Red
-- $1/4\lambda$ Slip
--Quartz Wedge

Polarizer, Sliding In/Out of Optical Path, On Top Of Collector



-- λ Slip First Class Red
-- $1/4\lambda$ Slip
--Quartz Wedge

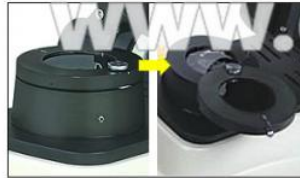

Polarizer, Sliding In/Out of Optical Path, On Top Of Collector



A15.1016 Polarizing Microscope

Rotatable Analyzer With Gradation 0° - 90°
-- λ Slip First Class Red
-- $1/4\lambda$ Slip Quartz Wedge

www.cnoec.com

Round Stage $\Phi 135\text{mm}$, 360° Rotatable



Opto-Edu (Beijing) Co., Ltd.

☎ 0086 13911110627

✉ sale@optoedu.com

🌐 cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China