



A43.1901 Digital Video Colposcope Universal Stand

Our Product Introduction

for more products please visit us on cnoec.com

Basic Information

- Place of Origin: China
- Brand Name: CNOEC, OPTO-EDU
- Certification: CE,
- Model Number: A43.1901
- Minimum Order Quantity: 1 pc
- Price: FOB \$1~1000, Depend on Order Quantity
- Packaging Details: Carton Packing, For Export Transportation
- Delivery Time: 5~20 Days
- Payment Terms: T/T, West Union, Paypal
- Supply Ability: 5000 pcs/ Month



Product Specification

- Observation Angle: 45° Tilted Binocular Observation Unit
- Eyepiece: 16x
- Interpupillary Distance: 50mm-80mm
- Diopter: $\pm 5D$
- Objective Focal Length: F=300mm (M58x0.75mm)
- Magnification Changer: 5-step Magnification Change: 0.4x, 0.66x, 1x, 1.5x, 2.5x
- Total Magnification: 2.7x, 4.4x, 6.6x, 10x, 16.6x
- Linear Field: 90mm, 55mm, 37mm, 24mm, 15mm
- Fine Focusing Distance: 20mm
- Floor Stand & Arm: Universal Stand,
- Highlight: **CE Digital Video Colposcope,
Universal Stand Digital Video Colposcope**

Product Description

A43.1901 Digital Video Colposcope Universal Stand

Table Of Specification For A43.1901 & A43.1902 Series Colposcope		
Item	A43.1901	A43.1902
Observation Angle	45° Tilted Binocular Observation Unit	
Eyepiece	16x	
Interpupillary Distance	50mm-80mm	
Diopter	±5D	
Objective Focal Length	F=300mm (M58x0.75mm)	
Magnification Changer	5-step Magnification Change: 0.4x, 0.66x, 1x, 1.5x, 2.5x	
Total Magnification	2.7x, 4.4x, 6.6x, 10x, 16.6x	
Linear Field	90mm, 55mm, 37mm, 24mm, 15mm	
Beam Splitter	Included, 50:50	
Video Camera Adapter	C-mount Included, For 1/3" Digital Camera	
Fine Focusing Distance	20mm	60mm
Illumination System	Coaxial Illumination, 10W LED, Brightness Adjustable, Illumination>30000lx	
Power Supply	Power Supply Input: AC110~240V, 50/60Hz,	
	With Power Supply Adapter Output DC12V3A, 36W,	
	Power Cable Length 3m With Magnetic Ring, Output Connector: Coaxial 5.5x2.1mm Plug	
Filter	Built-in Green And Blue Filters	
Floor Stand & Arm	Universal Stand,	Straight Vertical Stand,
	Balancing Arm With Counterweight & Lock	No Balancing Arm,
	Short Column Fitted On Five-star Base With Casters	Long Column Height Adjustable Fitted On The 5-jaws Round Base With Casters
Optional Accessories	Demonstrator, Camera (A59.4910, A59.4231), SLR Camera Adapter, Monitor, Objective (F=230mm & F=280mm)	




OPTO-EDU (BEIJING) CO., LTD.

OPTO-EDU

F-1501 Wanda Plaza, No.18 Shijingshan Road, Beijing 100043, China
Tel:+8610 88696020 Fax:+8610 88696085


A43.1901

Colposcope, Step Zoom 0.4x,0.6x,1x,1.6x,2.5x, Universal Stand

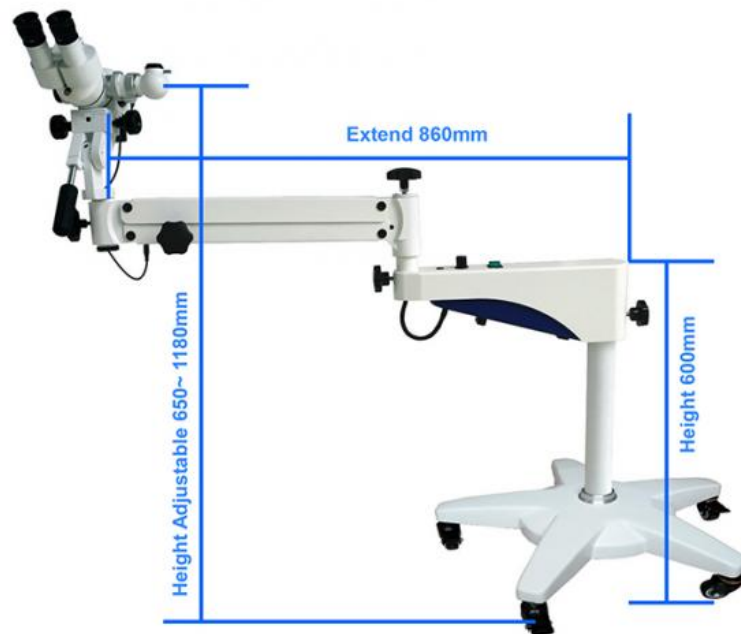

A43.1901 Optional Accessories

OPTO-EDU

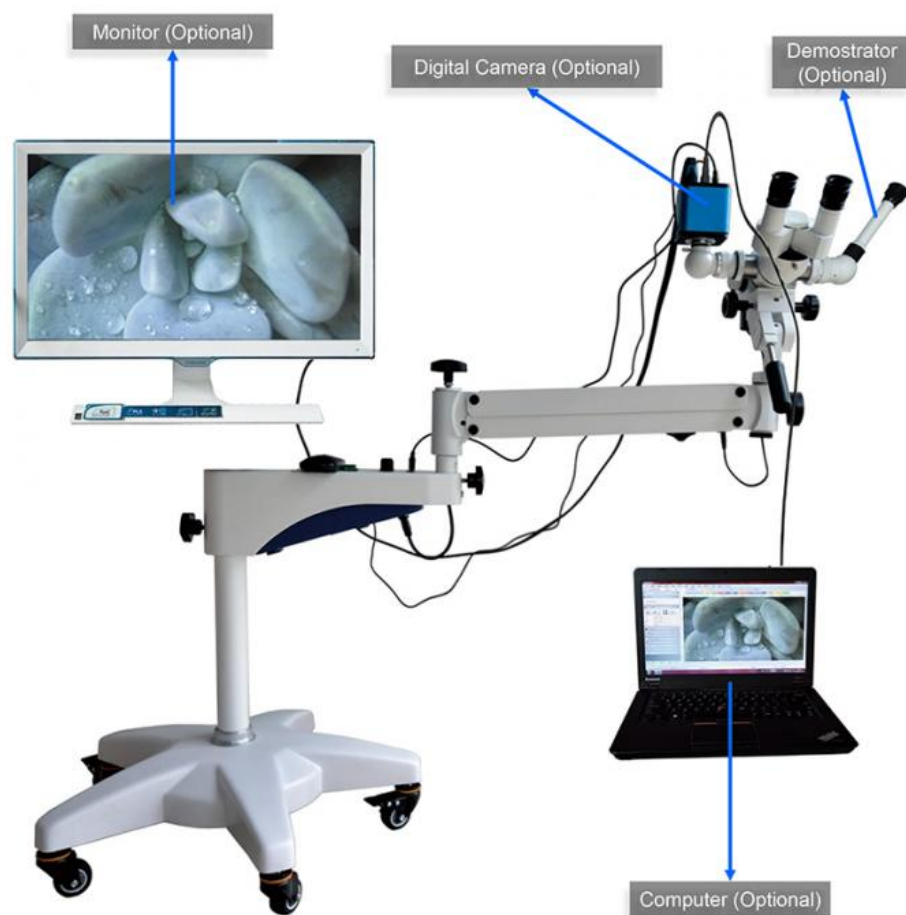



A43.1901



A43.1902

Table Of Optical Characteristics For A43.1901 & A43.1902 Series Colposcope (16X Eyepiece)

Focal Length Of Objective	F=300mm					F=230mm(Optional)					F=280mm(Optional)				
Magnification On Handwheel	2.5x	1.5x	1x	0.66x	0.4x	2.5x	1.5x	1x	0.66x	0.4x	2.5x	1.5x	1x	0.66x	0.4x
Total Magnification	16.6x	10x	6.6x	4.4x	2.7x	21.7x	13.2x	8.7x	5.7x	3.5x	17.8x	10.8x	7.1x	4.7x	2.9x
Linear Field (mm)	15	24	37	55	90	11	18	28	42	62	14	22	35	54	84
Resolution (LP/mm)	78	61	40	30	21	86	68	52	35	27	80	63	42	31	22

A43.1901 More Pictures

OPTO-EDU

Binocular Head 45° Tilted, IPD50-80mm, Diopter $\pm 5D$, Eyepiece WF16x



A43.1901 With Different Digital Camera (Optional) + Monitor/Computer (Optional)



A43.1901 With Different Digital Camera (Optional) + Monitor/Computer (Optional)







Opto-Edu (Beijing) Co., Ltd.

☎ 0086 13911110627

✉ sale@optoedu.com

🌐 cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China