



Wide Field Plan WF10X Biological Microscope A11.1316 With Double Layers Mechanical Stage

Our Product Introduction

for more products please visit us on cnoec.com

Basic Information

- Place of Origin: China
- Brand Name: CNOEC, OPTO-EDU
- Certification: CE, Rohs
- Model Number: A11.1316
- Minimum Order Quantity: 1 pc
- Price: Negotiation
- Packaging Details: Carton Packing, For Export Transportation
- Delivery Time: 5~20 Days
- Payment Terms: T/T, West Union, Paypal
- Supply Ability: 5000 pcs/ Month



Product Specification

- Eyepiece: Wide Field Plan WF10X
- Certification: CE|Rohs
- Focusing System: Coaxial Coarse And Fine Adjustment
- Illumination: LED
- Magnification: 40X-1000X
- Product Name: Biological Microscope
- Highlight: **high powered microscope ,
biological compound microscope**



More Images



Product Description

Wide Field Plan WF10X Biological Microscope A11.1316 With Double Layers Mechanical Stage

Introduction

A11.1316 is the best value entry level biological microscope, especially appreciated when budget is limited. All models are equipped with an integrated X-Y mechanical stage are supplied with an Abbe N.A. 1.25 condenser. Optional WF16X, WF20X eyepiece enable utmost 2000X optical magnification.



Specification

A11.1316 Biological Microscope		M	B	T
Head	Monocular Head, 30° Inclined, 360° Rotatable	√		
	Sliding Binocular Head, Interpupillary Distance: 55-75mm, 45° Inclined, 360° Rotatable, Left & Right Diopter Adjustable		√	
	Sliding Trinocular Head, Interpupillary Distance: 55-75mm, 45° Inclined, 360° Rotatable, Left & Right Diopter Adjustable			√
Eyepiece	Wide Field Plan WF10X	√	√	√
Objective	Achromatic 4X	√	√	√
	Achromatic 10X	√	√	√
	Achromatic 40X(S)	√	√	√
	Achromatic 100X(S,Oil)	√	√	√
Nosepiece	Quadruple Nosepiece	√	√	√
Stage	Double Layers Mechanical Stage 140X155mm/75X50mm	√	√	√
Condenser	Abbe Condenser N.A.1.25, Φ2-Φ30mm Iris Diaphragm, with Φ32 Filter Holder	√	√	√
Focusing	Coaxial Coarse & Fine Focusing Adjustable Mechanism 25mm, Precision 0.002mm	√	√	√
Light Source	LED Lamp 1W, Brightness Adjustable	√	√	√
Weight	Gross Weight 7Kg, Net Weight 6Kg	√	√	√
Package	1 Set/ Carton, 33cm*28cm*43cm, Wood Cabinet Is Optional	√	√	√

A11.1316 Optional Accessory		
Eyepiece	Wide Field Plan WF16X	A51.1303-16

PRODUCT DISPLAY



 A11.1316-B Monocular Head Biological Microscope



PRODUCT DETAILS



01 HEAD

Monocular, binocular, trinocular and teaching head available for choose. All heads inclined at 30°, rotatable 360°. For binocular and trinocular, both eyepiece tube diopter adjustable $\pm 5^\circ$.

02 NOSEPIECE

Quadruple nosepiece with precise click stop. Equipped with 45mm achromatic CF objective 4X, 10X, 40Xs, 100Xs (oil).




03 WORKING STAGE

Intergrated double layer mechanical stage, size 140x155mm, moving range 75x50mm.

04 FOCUSING

Coaxial coarse & fine focusing adjustment, with tension adjustment and Safety Stop Ring. Fine division 0.002mm, coarse stroke 36mm.

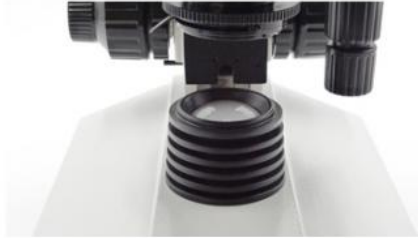

05 LIGHT SOURCE

High brightness 1W LED, brightness adjustable, with internal rechargeable batteries and external 100-240V battery charger/mains adapter.



Abbe Condenser N.A.1.25



Light Source Collector

RELATED PRODUCTS

Application

Biological microscope is widely used to observe biological slices, biological cells, bacteria and living tissue culture, fluid precipitation observation and research, while observing other transparent or translucent objects and powder, fine particles and other objects.



Opto-Edu (Beijing) Co., Ltd.



0086 13911110627



sale@optoedu.com



cnoec.com

F-1501 Wanda Plaza, No. 18 Shijingshan Road, Beijing 100043, China